

Cherrypicking with Source and Destination

 protocols.opentrons.com/protocol/7a3e4b

Author

Opentrons

This protocol utilizes a user-specified CSV to transfer sample using a variety of labware. With this protocol you have the ability to choose between the following:

- Starlab Deepwell Plate
- Eppendorf 2mL Tube in Tube Rack
- Eppendorf 1.5mL Tube in Tube Rack
- NEST 96-Well, 200µL Flat
- NEST 96-Well, 100µL PCR
- BioRad 96-Well, 200µL PCR
- Corning 96-Well, 360µL Flat
- Corning 384-Well, 112µL Flat
- USA Scientific 96-Deepwell, 2.4mL

With the destination labware providing the same options as source labware.

Explanation of complex parameters below:

Volumes CSV : Here, you should upload a .csv file formatted in the following way, being sure to include the header line:

Source Well | **Volume** | **Destination Well**


Note about CSV: The first column is dedicated to the header of the CSV; thus, the first row of data should occupy cells 'A2', 'B2', and 'C2' in the document.

- **Pipette Model** : Select which pipette you will use for this protocol.
- **Pipette Mount** : Specify which mount your single-channel pipette is on (left or right).
- **Source Labware Type** : Select which (destination) labware you will use for this protocol.
- **Destination Labware Type** : Select which (destination) labware you will use for this protocol.
- **Tip Use Strategy** : Select whether you would like to use the same tip or a new tip for each transfer.

Pipettes

Opentrons Single-Channel Pipette and corresponding Tips

Deck Setup

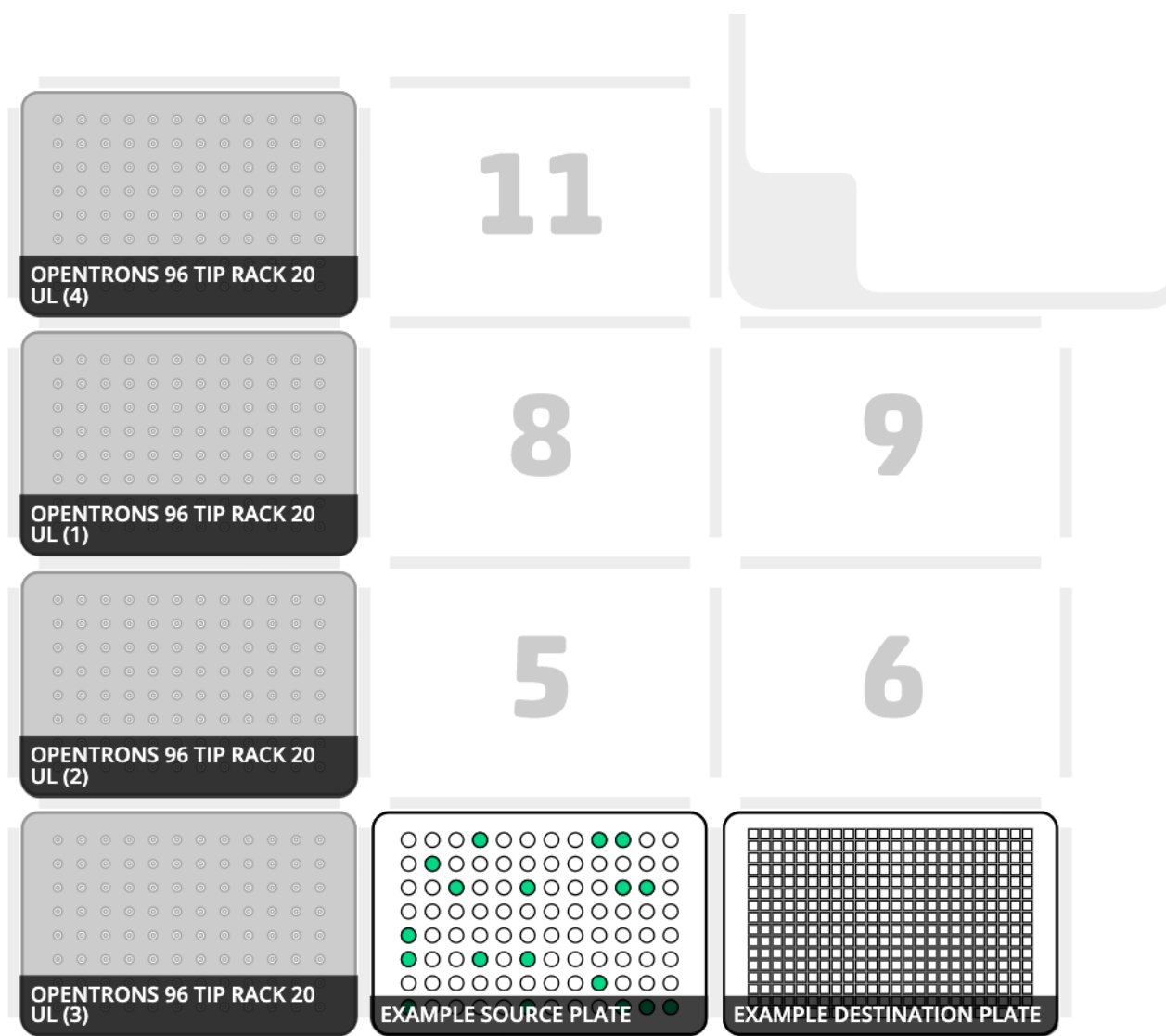

Protocol Steps

1. Pipette will aspirate a user-specified volume at the designated labware and well according to the imported csv file.
2. Pipette will dispense this volume into user-specified labware and well according to the imported csv file.
3. Steps 1 and 2 repeated over the duration of the CSV file.

OT-2

Download Protocol:

Process

1. Input your protocol parameters above.
2. Download your protocol and unzip if needed.

3. Upload your custom labware to the OT App by navigating to **More** > **Custom Labware** > **Add Labware** , and selecting your labware files (.json extensions) if needed.
4. Upload your protocol file (.py extension) to the OT App in the **Protocol** tab.
5. Set up your deck according to the deck map.
6. Calibrate your labware, tiprack and pipette using the OT App. For calibration tips, check out our [support articles](#).
7. Hit 'Run'.

Additional Notes

If you have any questions about this protocol, please contact the Protocol Development Team by filling out the [Troubleshooting Survey](#).